

eVO 350
eVO 450

OPHTHALMIC DIAGNOSIS



POWER OF PERFORMANCE AND PRECISION



With dedication to performance at an affordable price, these slit lamps were built for providing high definition imagery. Ergonomic configuration combined with superb optical systems and fluid motion, eVO350 and eVO450 provide reliable diagnosis. Both the configurations are integrated with highly efficient illumination system, an accurate apochromatic optical system with a platform fully expandable to adapt a host of accessories for extended performance.



mounted on robust and durable metal base plate, these systems are easy to install and configure to any office environment. these versatile slit lamps are ideal for efficient space management and optimum results.


sharing a common optical platform, eVO350 and eVO450 come built with two distinctly different illumination systems providing ideal slit lamp solutions for diverse needs. generating razor sharp images, superb contrast and high resolution, these systems bring to light the finest details of structures leading to accurate results.

EFFICIENCY AND FLEXIBILITY IN AN AFFORDABLE PACKAGE

the eVO 350 and eVO 450 introduce greater value for the office, thanks to easy adaptation of beam splitter with a built in camera port and focus adjust. this unique configuration of beam splitter keeps the magnichanger's footprint compact, ensuring comfort for the user and easy parfocalisation of camera with eye tubes.

a unique swing in Wratten filter makes the system most convenient for optometry. a goldmann tonometer (R-type and t-type) accessory with plug-and-play mounting promises precision in alignment and contact for accurate iOp readings. an optional diffuser filter kit can be introduced in the optical path for more enhanced region-of-interest illumination for brilliant photography.

all optics are lead-free and contact surfaces are treated with antimicrobial paint. patient chinrest is made of non-toxic and anti-irritation materials, promising safety for both professionals and patients alike.

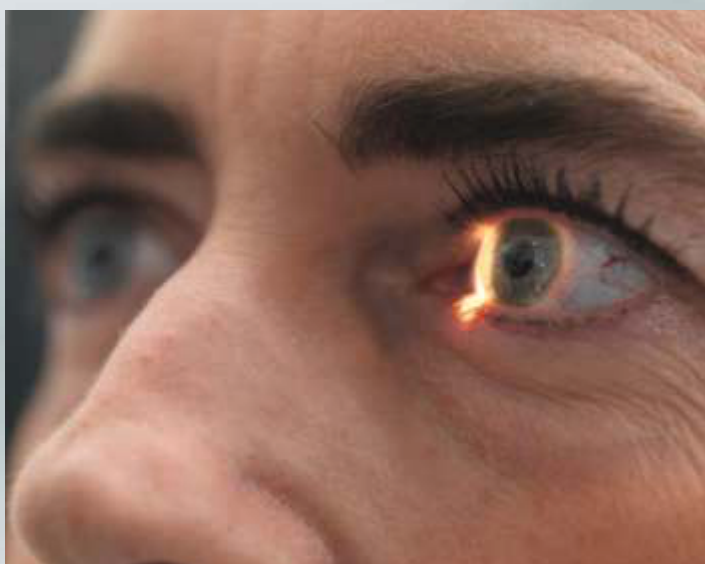


laBOMED maxlitetm coated optics yield crisp and detailed images that are true in fidelity and artifact-free. these proprietary coatings nearly eliminate internal reflections, while also making surfaces more reliably fog and scratch resistant.

eVO 350 / eVO 450

SPECIFICATIONS

Microscope	eVO 350	eVO 450
Magnification	Galilean	Galilean
Eyepieces	5 step: 6.5X, 10X, 16X, 25X, 40X	5 step: 6.5X, 10X, 16X, 25X, 40X
IPD range	WF 12.5X/18 mm, diopter adj. ± 5 mm	WF 12.5X/18 mm, diopter adj. ± 5 mm
Slit length/width	48.5 - 77 mm	48.5 - 77 mm
Slit apertures	0 - 13.5 mm	0 - 12 mm
Slit rotation	0.6, 5.8, 9, 13.5 mm and continuous	0.3, 1, 3, 5, 9, 12 mm and continuous
Slit tilt	1.5 to 12.0 mm	1.5 to 12.0 mm
Filters	0 - 180°	0 - 180°
Movement ranges:	-	20°
Longitudinal (in/out)	Red free, Cobalt blue, Heat absorbing	Red free, Cobalt blue, Heat absorbing,
Lateral (left/right)	Yellow filter (Wratten)	Neutral density, Yellow filter (Wratten)
Vertical (up/down)	100 mm	100 mm
Chin rest range	107 mm	107 mm
Illumination	30 mm	30 mm
Table base size	80 mm	80 mm
Weight	6V 20W Halogen	6V 20W Halogen
Optional accessories	465 mm X 316 mm	465 mm X 316 mm
	25 kg (packed)	24 kg (packed)
	Tonometer, Beam splitter, CCD, Diffuser,	Tonometer, Beam splitter, CCD, Diffuser,
	Motorized table	Motorized table

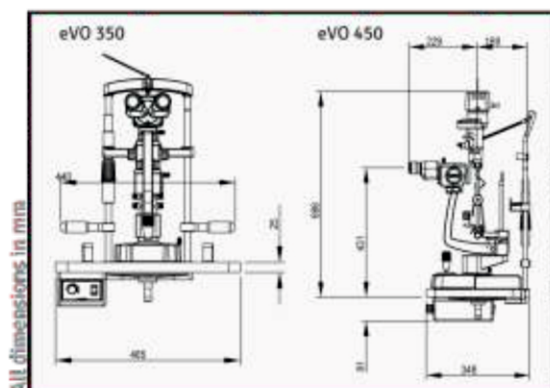

DESCRIPTION

eVO 350
eVO 450

CONFIGURATION

5 step Slit Lamp
5 step Slit Lamp

CATALOG NO.

8126300
8127300



LABOMED

Labo America, Inc.
920 Auburn Court
Fremont, CA 94538



Authorized Sales Partner:

Mex
ENTERPRISES

We Share in Healthcare...
ISO 13485 :2016 & ISO 9001 :2015 Certified Company

HEAD OFFICE

MEX ENTERPRISES
ISO 13485 :2016 & ISO 9001 :2015 Certified Company
Plot No. 50/IB-50/2B, Rajharsh
Colony, Near CI Square Shopping
Mall, Kolar Road, Bhopal, M.P
India-462042

DELHI SALES OFFICE

MEX ENTERPRISES
ISO 13485 :2016 & ISO 9001 :2015 Certified Company
Building Name - Classic's, T-1895/1,
Upper Ridge Road, Karol Bagh, New
Delhi India- 110005

RAIPUR SALES OFFICE

MEX ENTERPRISES
ISO 13485 :2016 & ISO 9001 :2015 Certified Company
House No.-10/806, ShivanandNagar,
Jhanda Chowk, Raipur Chhattisgarh
India-492008

Ph.: +91 755 4256391, +91 9302778762, +91 9111429222, Email: info@mediequipbazar.com, export@mediequipbazar.com

Website: www.mediequipbazar.com