


MEDIX

MIT-500

**MOTORIZED INSTRUMENT TABLE
WITH SLIDING TOP AND DRAWER**

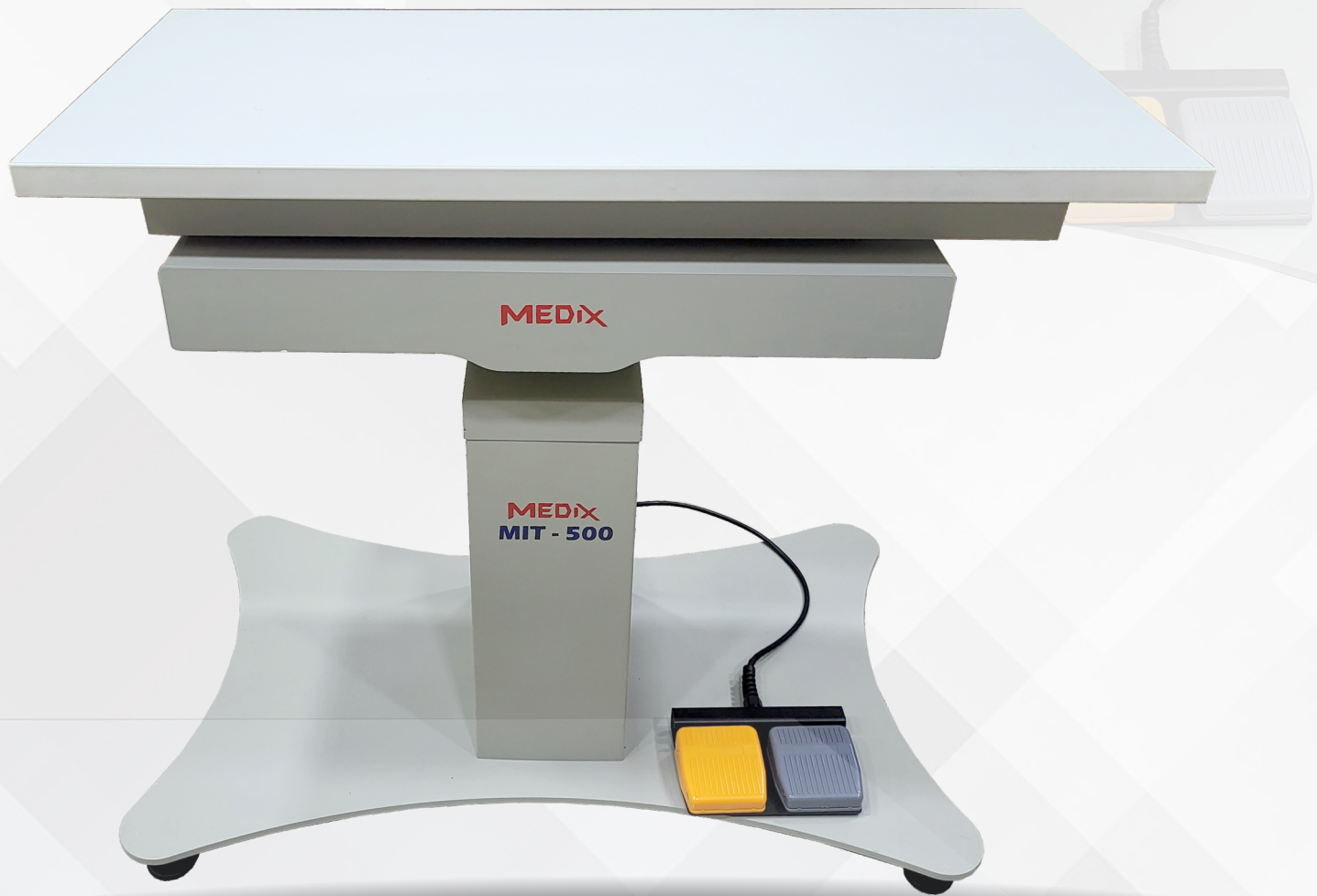


MEDIX


www.mediequipbazar.com

Motorized Instrument Table

- Very high quality, elegant and sturdy designed table for medical office, hospital, clinic and optical stores uses.
- Sliding top for two instruments makes it easy to operate with two different instruments without moving patient or operator.
- Motorized up-down movement with foot switch control provides hands free table height adjustment during procedures.
- High lifting capacity makes it suitable to use with heavy or large instruments as well.
- Imported silent DC motor with DC power supply makes the operation noiseless and safe.
- Built-in large drawer provides adequate space for trial set.

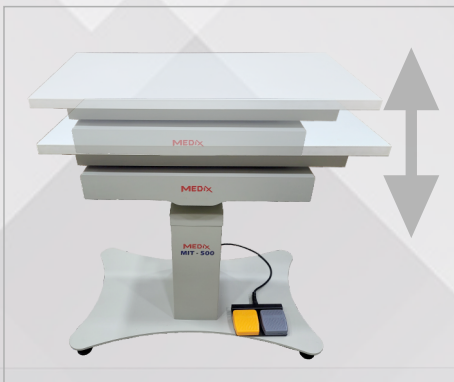

Motorized Up-Down

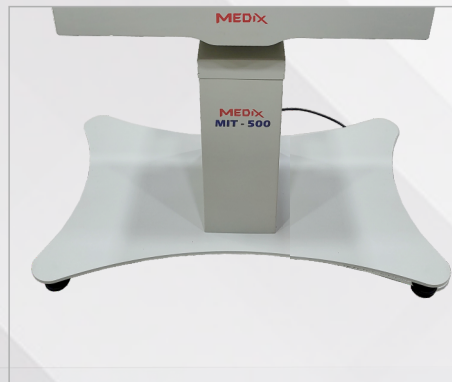


Stable Base



Sliding Switch Top



Double Instrument Top

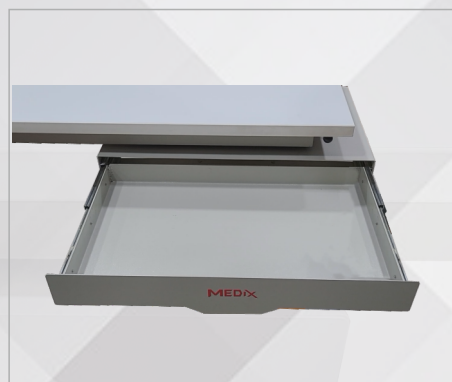

Large Drawer



Easy Foot Switch Control



TECHNICAL SPECIFICATION

Maximum Height:	900 mm
Minimum Height:	720 mm
Height Adjustment Control	Foot Switch Control
Vertical Adjustment Range:	180 mm
Total Length:	800 mm
Total Width:	500 mm
Table Top Sliding Range:	400 mm
Sliding Mechanism Control:	Locking Lever pull and release
Table Top Length:	800 mm
Table Top Width:	400 mm
Drawer Dimensions:	540 mm (L) x 340 (W) mm
Power:	100~240 AC 50Hz
Auxiliary Connection:	1 Auxiliary Socket
Weight Capacity:	75 Kilograms
Gross Weight:	70 Kilograms
Standard Accessories:	Foot switch

MEDIx



+91 9302 778 762
+91 9111 429 222



export@mediequipbazar.com
info@mediequipbazar.com



+91 9302 778 762
+91 9111 429 222

REGISTERED OFFICE

MEX ENTERPRISES

ISO 13485 :2016 & ISO 9001 :2015 Certified Company

Plot No. 50/1B-50/2B, Rajharsh Colony,
Near CI Square Shopping Mall, Kolar Road,
Bhopal, M.P.-462042

SALES OFFICE

MEX ENTERPRISES

ISO 13485 :2016 & ISO 9001 :2015 Certified Company

Building Name - Classic's, T-1895/1,
Upper Ridge Road, Karol Bagh,
New Delhi - 110005

SALES OFFICE

MEX ENTERPRISES

ISO 13485 :2016 & ISO 9001 :2015 Certified Company

House No.-10/806, Shivanand Nagar,
Jhanda Chowk, Raipur CG-492008



MANUFACTURING UNIT

Plot No. 23-B, Badiyakhedi
Industrial Park Phase I,
Badiyakhedi, Sehore, M.P.,
India - 466 001



MANUFACTURING UNIT

Plot No.-40, KI Ansal
Green, Kolar Road
Bhopal-462042



MANUFACTURING UNIT

Plot No. - 79, Commercial
Shade, Near D-Mart,
Bhopal, 462042



www.mediequipbazar.com

MEDIx



www.mediequipbazar.com